

Earmould insertion - step by step guide




Hearing aid

Earmould

Red = right aid

Blue = left aid

Step 1 - hold the edge of the earmould between the index finger and thumb

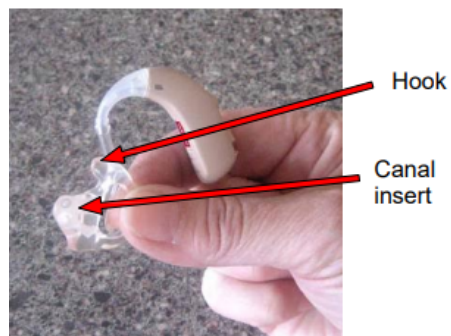

Step 2 – lift the earmould to the ear using the other hand to guide




Step 3 – slide backwards & twist the mould back to fit hook under top flap



Step 4 – gently push the mould into the ear & check the hook is under top flap



Step 5 – lift the hearing aid over the top of the ear to keep it in place

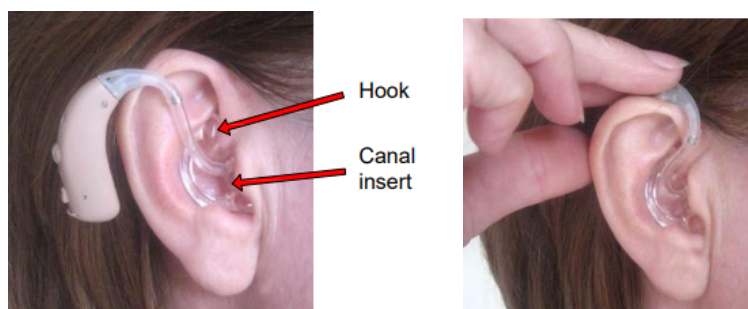

Wipe the earmould with a damp cloth or wet wipe every night after use. Please refer to the Hearing Aid Maintenance advice sheet for information on cleaning and checking the hearing aid and mould alongside retubing if required.